



Beijing Emarbio Science & Technology Co.,LTD

Anti-HSP90 α Rabbit Polyclonal Antibody

Catalog No.	EM42151-01	EM42151-02
Size	50 μ l	100 μ l
Species	Human Mouse Rat	<i>Molecular</i> 95KD

Background: Hsp90 (heat shock protein 90) is a chaperone protein that assists other proteins to fold properly, stabilizes proteins against heat stress, and aids in protein degradation. In mammalian cells, there are two or more genes encoding cytosolic Hsp90 homologues, with the human Hsp90 α showing 85% sequence identity to Hsp90 β .

Specificity: Antibody can detect endogenous HSP90 α protein.


Application: WB 1:1000-2000

Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid, 1mg/ml

Storage buffer: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Storage: Store at -20°C. Do not aliquot the antibody.

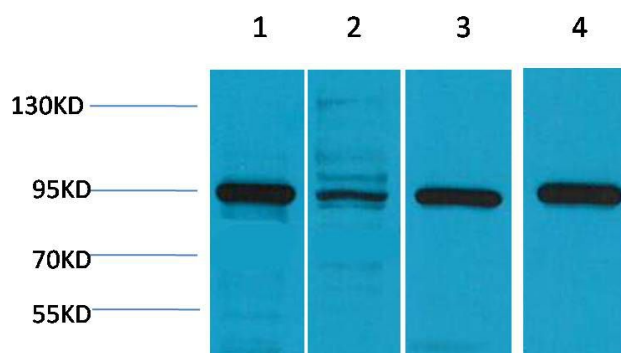

Western blot analysis of 1) HeLa, 2) 3T3, 3) Mouse Brain, 4) Rat Brain, with HSP90 α Rabbit pAb diluted at 1:2,000.