



Anti-cMyc Tag Rabbit Polyclonal Antibody

Catalog No.	EM34014-01	EM34014-02
Size	50 μ l	100 μ l

c-Myc-tag antibody is part of the Tag series of antibodies, the best quality in the research. Myc protein is a transcription factor that activates expression of a great number of genes through binding on consensus sequences (Enhancer Box sequences (E-boxes)) and recruiting histone acetyltransferases (HATs). A recent study demonstrated that temporary inhibition of Myc selectively kills mouse lung cancer cells, making it a potential cancer drug target.

Molecular Weight: N/A


Application: WB 1:5,000 IP 1:200 Optimal dilutions should be determined by the end user.

Specificity: The C-Myc tag antibody can detect C-Myc fusion proteins.

Form: Liquid ,1mg/ml

Storage buffer: PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol

Storage: Store at -20°C. Do not aliquot the antibody.

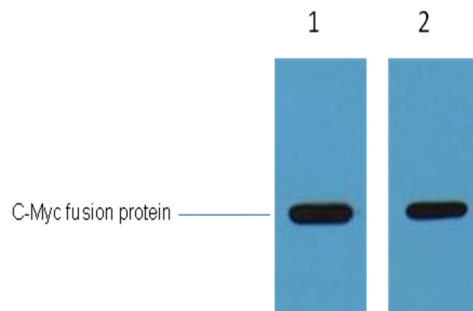

1 μ g C-Myc fusion protein+ Primary antibody dilution at
1、 1:5,000 2、 1:10,000

