



anti-Histone H3 Rabbit Polyclonal antibody

Catalog No.	EM12348-01	EM12348-02
Size	50µl	100µl

Species reactivity: **Human Mouse Rat** Molecular Weight: 17KD

Background: H3 Core component of nucleosome. Nucleosomes wrap and compact DNA into chromatin, limiting DNA accessibility to the cellular machineries which require DNA as a template. Histones thereby play a central role in transcription regulation, DNA repair, DNA replication and chromosomal stability. DNA accessibility is regulated via a complex set of post-translational modifications of histones, also called histone code, and nucleosome remodeling. The nucleosome is a histone octamer containing two molecules each of H2A, H2B, H3 and H4 assembled in one H3-H4 heterotetramer and two H2A-H2B heterodimers.

All Names: H3 histone family, member A; H3/A; H31; H3FA; H3FB; H3FC; H3FD; H3FF; H3FH; H3FI; H3FJ; H3FK; H3FL; HIST1H3A; HIST1H3B; HIST1H3C; HIST1H3D; HIST1H3E; HIST1H3F; HIST1H3G; HIST1H3H; HIST1H3I; HIST1H3J; histone 1, H3a; histone cluster 1, H3a; Histone H3.1;

Immunogen: A synthesized peptide derived from human Histone H3.

Purification: Immunogen affinity purified


Application: WB 1:500-1:3000 IHC: 1:50-1:200 IF: 1:100-1:500

Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid ,1mg/ml

Storage buffer: Rabbit IgG in phosphate buffered saline , pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol. Store at -20 ° C. Stable for 12 months from date of receipt

Storage: Store at -20°C. Do not aliquot the antibody.

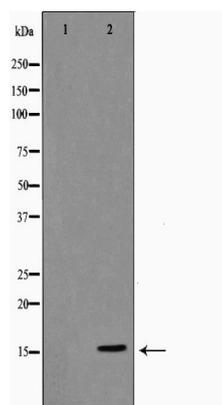

Western blot analysis on COLO205 cell lysate using Histone H3 antibody, The lane on the left is treated with the antigen-specific peptide.