



EMAR Beijing Emarbio Science & Technology Co.,LTD

Anti-Cytokeratin 6 Mouse Monoclonal Antibody(5B3)

Catalog No.	EM42089-01	EM42089-02
Size	50 μ l	100 μ l
Species	Human	Molecular Weight: 138KD

Background: Cytokeratins are proteins of keratin-containing intermediate filaments found in the intracytoplasmic cytoskeleton of epithelial tissue.

Specificity: Cytokeratin 6 Mouse Monoclonal antibody detects endogenous Cytokeratin 6 proteins.


Application: WB: 1:1,000-3000 IHC:1:200

Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid ,1mg/ml

Storage buffer: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Storage: Store at -20°C. Do not aliquot the antibody.

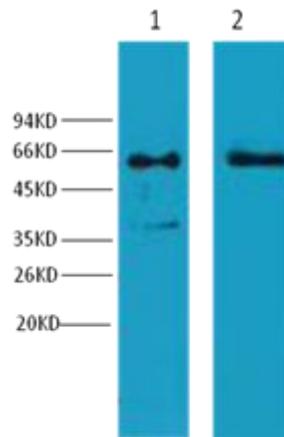

Western blot analysis of 1) HeLa, 2) HepG2 ,with Cytokeratin Mouse mAb diluted at 1:2,000.

EMAR

EMAR