



anti-Bax(Ab-167) Rabbit Polyclonal antibody

Catalog No.	EM51052-01	EM51052-02
Size	50µl	100µl
Species reactivity:	Human Mouse Rat	Molecular Weight:

All Names:Apoptosis regulator BAX, membrane isoform alpha, BAXA, Bcl2-associated X protein

Immunogen:The antiserum was produced against synthesized non-phosphopeptide derived from human Bax around the phosphorylation site of threonine 167 (F-G-T^P-P-T).

Purification:The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.


Application: WB 1:500-1:3000 IF 1:100-1:500 ELISA 1:1000

Optimal dilutions/concentrations should be determined by the end user.


Form: Liquid ,1mg/ml

Storage buffer: Rabbit IgG in phosphate buffered saline , pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.Store at -20 ° C.Stable for 12 months from date of receipt

Storage: Store at -20°C. Do not aliquot the antibody.

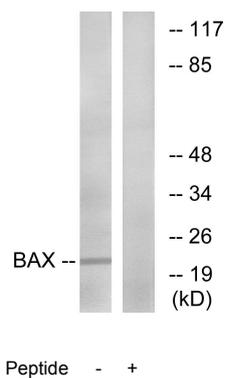

Western blot analysis of extracts from HepG2 cells, using Bax (Ab-167) antibody.

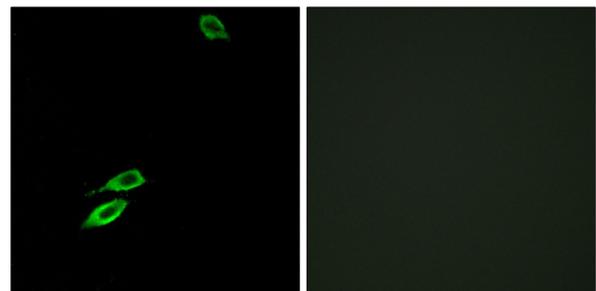

Immunofluorescence analysis of LOVO cells, using Bax (Ab-167) antibody.