



EMAR Beijing Emarbio Science & Technology Co.,LTD

Anti-FOXO1A Rabbit Polyclonal Antibody

Catalog No.	EM52105-01	EM52105-02
Size	50 μ l	100 μ l
Species reactivity:	Human Mouse Rat	Molecular Weight: 70KD

Background: This gene belongs to the forkhead family of transcription factors which are characterized by a distinct forkhead domain. The specific function of this gene has not yet been determined; however, it may play a role in myogenic growth and differentiation.

All Names: FKHR, FOXO1, Forkhead box protein O1A, Forkhead in rhabdomyosarcoma

Immunogen: A synthesized peptide derived from human FOXO1A

Specificity: FOXO1A Antibody detects endogenous levels of total FOXO1A


Application: WB: 1:1,000-3000 IHC 1:50-1:200 IF/ICC:1:100-1:500

Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid ,1mg/ml

Storage buffer: PBS, pH 7.4, containing 0.02% *sodium azide* as Preservative and 50% Glycerol.

Storage: Store at -20°C. Do not aliquot the antibody.

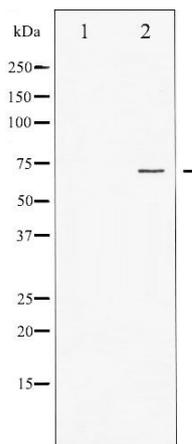

Western blot analysis of FOXO1A expression in Serum treated HeLa whole cell lysates, The lane on the left is treated with the antigen-specific peptide.

EMAR

EMAR